



Designed by Scientists for Scientists

Scientific Behavioral Testing Systems

SAN DIEGO INSTRUMENTS

Product Catalog



The World's Most Powerful, Flexible and Easy to Use Behavioral Testing

Welcome



San Diego Instruments (SDI) has served the scientific community for over 30 years as a comprehensive resource for the design, manufacture and distribution of scientific testing systems used in human and animal behavioral studies. Utilized in laboratories and cited in research papers worldwide, SDI systems represent the industry benchmark for quality and longevity. Our commitment to developing quality products that stand the test of time is matched only by our dedication to excellent customer service. We take pride in our ever-growing core of loyal clientele.

We are proud to be poised on the cutting-edge of industry trends. Our continual investment in design and development keeps our products one step ahead of the rapidly evolving needs of our clientele. SDI behavioral testing systems afford you the utmost in quality and performance, giving you the edge in an industry where Power, Flexibility and Ease of Use are everything.

Power

The power of SDI test systems translates directly into speed and efficiency in the laboratory. Taking full advantage of modern computers, our data organization and management software provides extensive test options and the ability to configure multiple test stations, allowing rapid testing of a large number of subjects. One-touch exporting of data to statistical packages simplifies the process by eliminating the hassle of cutting and pasting.

Flexibility

Specifically designed with flexibility in mind, SDI systems are easily customized to your unique testing parameters. System configurations can be easily expanded over time. Intuitive yet powerful features cater to users from the educated non-specialist to the sophisticated behaviorist. Our systems support multiple test paradigms without requiring costly add-on kits or additional software. What's more, several of our systems are capable of running on a desktop or laptop computer. (Soon all SDI systems will offer this ability).

Ease of Use

All the power and flexibility in the world, however, is nothing without ease of use. All SDI test systems are self-installable in less than one hour. Simple, menu-driven steps allow you to set up and run tests without any programming. Our consistent point-and-click user interfaces minimize technician training and none of our systems require pre-training of the subject.

SDI behavioral testing systems are available within and outside of the United States. Internationally, our products are sold through a network of distributors. If you have any questions or would like to request a quote please call (858) 530-2600 or email us at sales@sandiegoinstruments.com. Visit us online at www.sandiegoinstruments.com.

SDI would like to thank all the scientists worldwide who have contributed to the success of our behavioral testing systems – our systems are truly Designed by Scientists for Scientists. We look forward to assisting you in the advancement of your research.

Table of Contents

STARTLE RESPONSE SYSTEMS

SR-LAB™	4
---------	---

MOTOR & SENSORY SYSTEMS

ROTOR-ROD™	6
Animal Grip Strength	7

LEARNING SYSTEMS

GEMINI™	7
Freeze Monitor™	9

ACTIVITY SYSTEMS

Photobeam Activity System - Home Cage	11
Photobeam Activity System - Open Field	11
Rotometer	12

PLACE PREFERENCE

TREMOR MONITOR™

VIDEO TRACKING SYSTEMS

ANY-Maze™	14
-----------	----

MAZES & ENCLOSURES

Barnes Maze, Y-Maze, Elevated Plus Maze, Radial Arm Maze, T-Maze, Zero Maze, Multiple Unit Open Field, Single Unit Open Field, and Water Maze.

HUMAN TESTING SYSTEMS

SR-HLAB™	19
Eyeblink Conditioning System	20
ERP-LAB™	21
Human Testing Combination Systems	21

TEST ACCESSORIES


AIRSTIM™	22
AIRSTIM™ with Software	22
Portable Air Puff Unit	22
Portable Air Puff Unit with Software	22
Animal Stimulator	22
Animal Shocker	22
Shock Level Tester	22
White Noise Generator	22

SDI CONFIGURED COMPUTERS

SDI GENERAL INFO

STARTLE RESPONSE SYSTEM

SR-LAB™

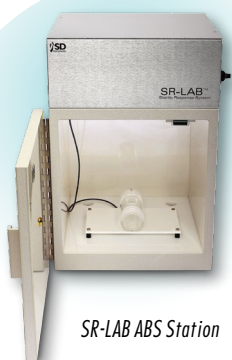

SR-LAB ABS Station

The SR-LAB™ Startle Response System – the world's most widely used system for startle reflex measurement and by far the most successful for fear potentiated startle and pre-pulse inhibition testing. SR-LAB provides a complete hardware and software solution for a wide variety of startle applications. Intuitive yet powerful features cater to users from the educated non-specialist to the sophisticated behaviorist.

SR-LAB takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. SR-LAB software controls virtually any combination of tones, noise bursts, lights, air puffs, background noise and foot shock options. Flexible data organization and management software provides the ability to configure multiple test stations, allowing rapid testing of a large number of subjects. Importantly, SR-LAB supports multiple test paradigms without requiring costly add-on kits or additional software.

• Maximum Stations: 16

For more detailed SR-LAB info:
www.sandiegoinstruments.com/srvideo


2325-0400	SR-LAB ABS SYSTEM Includes: 1 ABS cabinet, 1 animal enclosure, software with Analysis data analyzer, connection chassis for 16 stations with USB Interface and all cables & connectors.
2325-0401	SR-LAB ABS System - Additional Test Station
2325-0402	SR-LAB PRIME SYSTEM Includes: 1 Prime cabinet, 1 animal enclosure, software with Analysis data analyzer, connection chassis for 16 stations with USB Interface and all cables & connectors.
2325-0403	SR-LAB Prime System - Additional Test Station
See pg.14	DUAL USAGE SR-LAB/TREMOR MONITOR Includes: 1 ABS cabinet, 1 Tremor Monitor animal enclosure, 1 SR-LAB animal enclosure, Tremor Monitor software, SR-LAB software, connection chassis with USB Interface and all cables & connectors.

SR-LAB OPTIONAL ACCESSORIES

2325-0078	Standardization Unit
Additional Animal Enclosures:	
7500-0320	Small with UltraSensor Plate (3.5" L x 1.1" I.D.) (up to 25 G.)
7500-0321	Small/Medium with UltraSensor Plate (5" L x 1.5" I.D.) (up to 50 G.)
7500-0322	Medium (6" L x 2.25" I.D.) (Up to 150 G.)
7500-0323	Large (8" L x 3.5" I.D.) (Up to 450 G.) - standard rat
7500-0324	Extra Large (10" L x 5" I.D.) (Up to 800 G)
7500-0325	Small/Medium Animal Enclosure Slotted for Tethered Subjects
7500-0326	Medium Animal Enclosure Slotted for Tethered Subjects
7500-0327	Large Animal Enclosure Slotted for Tethered Subjects
7500-0328	Extra Large Animal Enclosure Slotted for Tethered Subjects



SR-LAB Enclosure, Base and Standardization Unit

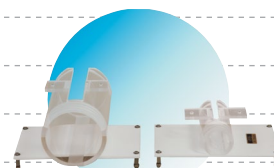


Slotted Animal Enclosures

SR-LAB STIMULI ACCESSORIES

Air:


- 2325-0368 Air Puff Kit - Per test station
- 2200-0101 Air Puff Distribution Box - Digital Port Connection
- 7500-0292 Air Puff Air Manifold(Required for multiple stations. Max: 8 Stations)
- 7500-0308 Air Puff Delivery Assembly(Air Nozzle, Animal Enclosure Mount and Tubing)

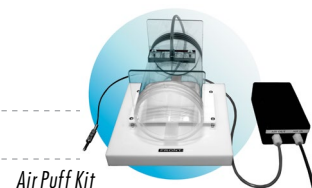

Air Puff Kit

Light:

- 7500-0210 PPI-LED Light Unit (Per Station)
- 7500-0243 High Intensity Light Bar(Per Station)
- 2200-0102 Lights Distribution Box - Digital Port Connection

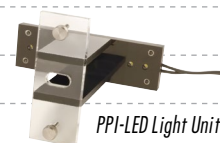


PPI-LED Light Unit

Acoustic Stimulus - Tone

- 7500-0350 Pure Tone Auditory System (software, amplifier for up to 4 test stations, and all cables and connectors)
- 7500-0329 Pure Tone Audio Amplifier (supports 4 additional test stations)
- 7500-0351 Pure Tone Speaker Assembly (Per Station)

Shock:

- 2325-0404 Potentiated Startle Kit
Includes: Shocker, Shock Grid, LED Signal Light, and all cables & connectors..
- 7500-0007 SR-LAB Small Shock Grid.
- 7500-0006 SR-LAB Small/Medium Shock Grid
- 7500-0005 SR-LAB Medium Shock Grid
- 7500-0004 SR-LAB Large Shock Grid
- 7500-0011 SR-LAB X-Large Shock Grid

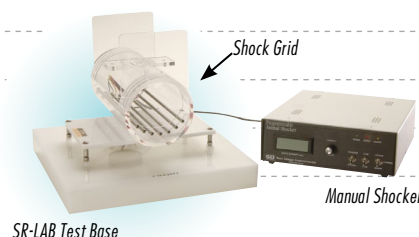

SR-LAB Test Base

Manual Shocker

SYSTEM AND SOFTWARE UPGRADES

- 6300-0000 SR-LAB Windows® Software Version Upgrade
- 2325-0380 SR-LAB DOS/PCI to SR-LAB Windows USB Upgrade



Animal Enclosures

MOTOR & SENSORY SYSTEMS

ROTOR-ROD™

ROTOR-ROD™
Animal Grip Strength



ROTOR-ROD Station with Rat Cart

The ROTOR-ROD™ System measures motor function, motor learning, coordination, and equilibrium in both rats and mice. ROTOR-ROD takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. ROTOR-ROD is designed with adjustable falling heights from 18-48 inches to employ 'fear of falling' as a natural motivator.

ROTOR-ROD is available in three configurations: a model for mice, a model for rats and a combination model for both rats and mice. Each model features four independently timed lanes that automatically record fall latencies and distances traveled. Secured enclosures prevent the subjects from escaping after a safe landing.. The user defines the segments for rod speeds; constant, accelerating or decelerating. The Rotor Rod software allows as many segments as will fit within the maximum time for a test.

For more detailed Rotor-Rod info:
www.sandiegoinstruments.com/rotorrodvideo

ROTOR-ROD SYSTEM

Includes: Software, USB to serial adapter and all cables & connectors.

2325-0405	ROTOR-ROD System - Mouse
2325-0406	ROTOR-ROD System - Rat. Includes Rat Cart
2325-0407	ROTOR-ROD System - Rat & Mouse. Includes: Rat Cart, Rat Rod and Mouse Rod
7500-0330	Rotor Rod System - Mouse Model to Mouse and Rat Model
7500-0331	Rotor Rod System - Rat Model to Mouse and Rat Model

ROTOR-ROD 8" DISKS UPGRADE KITS

2325-0256	Rotor Rod System - Upgrade Mouse to Rat and Mouse
2325-0258	Rotor Rod System - Upgrade Rat to Rat and Mouse

ANIMAL GRIP STRENGTH SYSTEM



*Grip Strength Platform with 2 gauges
for High Throughput*

SDI's Animal Grip Strength System measures hindlimb and forelimb grip strength. Gauge mountings are horizontally and vertically adjustable to accommodate both rats and mice. Two sets of walls for the center channel are included: tall for rats and short for mice. The white acrylic platform is easy to disassemble for cleaning. Models are available that include one or two gauges with PC interface for automatic recording. The Grip Strength Platform with two gauges is ideal for high throughput testing, allowing the user to test both hindlimb and forelimb strength in a single pass. The adjustable Grip allows the user to position the angle of the grip between 0 and 45 degrees to obtain the optimal results.

ANIMAL GRIP STRENGTH SYSTEM

All gauges come with a PC interface and all systems include software. The gauges can be used in manual mode without connecting to a computer if desired.

2325-0063	Platform with 1 Gauge and 1 Grip (Specify rat or mouse grip)
2325-0064	Platform with 2 Gauges and 2 Grips (Specify rat or mouse grips)

LEARNING SYSTEMS

GEMINI™

Freeze Monitor™

GEMINI™



GEMINI Test Station with Shock

For more detailed GEMINI info:
www.sandiegoinstruments.com/geminivideo


The GEMINI™ Avoidance System is designed for active and passive avoidance and learned helplessness testing in both rats and mice. GEMINI is available in two configurations: Shock or Shock and Airstim. GEMINI takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. Conversion from rats to mice is possible in less than five minutes. GEMINI provides maximum flexibility in test paradigms by utilizing features such as non-heating LED house lights, laboratory standard cue lights, user selectable auditory stimuli, infrared photobeams and a quiet auto door. Each enclosure starts and runs independently making continuous testing quick and easy. Importantly, GEMINI can be easily expanded to accommodate multiple test stations, allowing rapid testing of a large number of subjects. Supported avoidance paradigms are: Passive Avoidance (Classic, Trials to Criterion, Continuous), Active Avoidance (Trials to Criterion, Continuous) and Learned Helplessness.

• **Maximum Stations: 8**

2325-0301	GEMINI AVOIDANCE SYSTEM WITH SHOCK™ (Specify for rats or mice) Includes: Software, auto door, grid floor, LED house lights, cue lights, auditory stimuli and all cables & connectors.
2325-0302	GEMINI Avoidance System - Additional Test Station with Shock (Specify for rats or mice)
2325-0303	GEMINI AVOIDANCE SYSTEM WITH SHOCK AND AIRSTIM™ (Specify for rats or mice) Includes: Software, auto door, grid floor, LED house lights, cue lights, auditory stimuli and all cables & connectors.
2325-0304	GEMINI Avoidance System - Additional Test Station with Shock and AIRSTIM™ (Specify for rats or mice)
GEMINI SYSTEM OPTIONAL COMPONENTS & ACCESSORIES	
7500-0009	Auto Door for Passive Avoidance
7001-0127	Mouse "Start Box" for Passive Avoidance
7500-0018	Gemini Grid Floor for Rats
7500-0019	GEMINI Grid Floor for Mice
2325-0104	GEMINI Isolation Cabinet
2325-0361	GEMINI Avoidance System-Shock Upgrade PCI Model to USB Model. <i>See note #1</i>
2325-0363	GEMINI Avoidance System-Shock/Air Upgrade PCI Model to USB Model. <i>See note #1</i>
6300-0200	GEMINI Windows® Software Version Upgrade
2325-0358	Shock Level Tester

Note #1: Some existing systems can be upgraded to the USB version. Serial Numbers: GEM-233 and higher can be upgraded. This requires the system be returned to SDI for upgrading. The price assumes reusability of all electronic components, the addition of the required USB components, rewiring and GEMINI USB Software. If components cannot be reused the customer will be contacted with a revised quotation to do the upgrade. USB Version additional features: Sound Cue, White Noise or Specific Frequency (Variable between 4 and 19 KHz), LED House Light(no heat), Delayed Shock in Passive Avoidance and inclusion of Learned Helplessness paradigm.

FREEZE MONITOR SYSTEM™

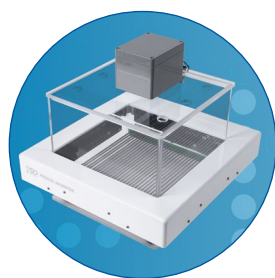

The Freeze Monitor™ System is designed to study cued and contextual fear conditioning. The subject can be forced into a freezing defensive posture by utilizing auditory cues, light cues and a shock grid. Freeze Monitor software provides the user with the ability to control the cue type, cue duration and scoring frequency. The 16 x 16 photobeam array precisely records the start and end of freezing episodes. Freeze Monitor software accurately reports the total number of freezing episodes and total time freezing for each session. Freeze Monitor takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use.

• **Maximum Stations: 4**

2325-0185	FREEZE MONITOR SYSTEM FOR RATS Includes: Software, acrylic test station, grid floor, light/sound stimulus unit and all cables & connectors.
2325-0188	Additional Freeze Monitor Test Station (Rats)

2325-0241 **FREEZE MONITOR SYSTEM FOR MICE**
Includes: Test station, grid floor, software, light/sound stimulus unit and all cables & connectors.

2325-0242 Additional Freeze Monitor Test Station (Mice)

FREEZE MONITOR OPTIONAL COMPONENTS & ACCESSORIES


7500-0195 Freeze Monitor Grid Floor (Mice)

7500-0196 Freeze Monitor Grid Floor (Rats)

7500-0176 Isolation Cabinet with Fan and Light


7500-0278 Freeze Monitor Context Box

2325-0358 Shock Level Tester

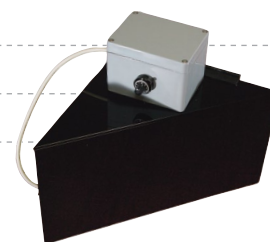

Freeze Monitor Context Box

PHOTOBREAM ACTIVITY SYSTEM - HOME CAGE

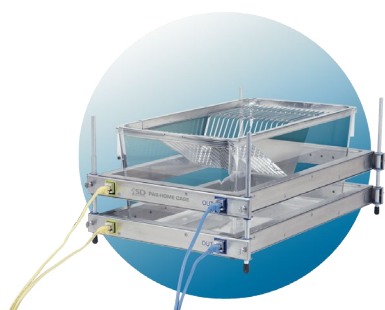

PAS - Home Cage with Rearing Option

For more detailed PAS Home Cage info:
www.sandiegoinstruments.com/pasvideo

SDI's Photobeam Activity System – Home Cage (PAS-HC) provides an affordable and flexible solution for activity testing. PAS-HC utilizes a laboratory's standard transparent cages, thus reducing the cost of a system. PAS-HC supports the rapid exchange of subjects in short studies, while affording the ability to run 24 hours a day for a full week without user intervention. PAS-HC utilizes a 4 x 8 photobeam configuration to track the subject's path in its home cage. Optional features include rearing frames that are infinitely adjustable to accommodate subjects of different sizes. PAS-HC takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. All study data is stored in a database for easy exporting to statistical packages. PAS-HC accurately records ambulation movements, fine movements, rearing (if option is present), time stamped (x,y) positions and central vs. peripheral counts. Through the use of the PAS Reporter utility, zone entries, time in zone, distance traveled, speed and resting time can also be reported. Powerful software features provide the ability to graphically view the subject's path.

• **Maximum Stations: 32 with or without rearing frames** (For older PCI interface PAS Home Cage Systems please contact San Diego Instruments to determine if your system can be supported, enlarged and/or converted to a USB interface.)

ACTIVITY SYSTEMS

Photobeam Activity System - Home Cage

Photobeam Activity System - Open Field

Rotometer

2325-0408 **PHOTOBREAM ACTIVITY SYSTEM - HOME CAGE**
Includes: Control Unit with USB interface, software, 4 x 8 mounting photobeam frame and all cables & connectors.

7500-0332 Photobeam Activity System - Home Cage Additional Test Station

PHOTOBREAM ACTIVITY SYSTEM - HOME CAGE OPTIONAL COMPONENTS & ACCESSORIES

7500-0333	Photobeam Activity System - Home Cage Rearing Frame
6300-0100	Photobeam Activity System - Home Cage Windows® Software Version Upgrade
6300-0107	Photobeam Activity System - Home Cage FDA Part 11 Compliance Software
6603-0100	Photobeam Activity System - Home Cage GLP Validation Document Set

PHOTOBREAM ACTIVITY SYSTEM - OPEN FIELD



PAS-Open Field with Rearing Option

For more detailed PAS Open Field info:
www.sandiegoinstruments.com/pasvideo


SDI's Photobeam Activity System – Open Field (PAS-OF) provides a powerful solution for applications in which qualitative analysis of the subject's locomotor patterns is necessary. PAS-OF utilizes a 16 x 16 photobeam configuration to track the subject's path within an acrylic enclosure. Optional features include rearing frames that are infinitely adjustable to accommodate subjects of different sizes and floor pokes. PAS-OF takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. PAS-OF accurately records and reports ambulation movements, fine movements, rearing (if option is present), pokes (if option is present), time stamped (x,y) positions and central vs. peripheral counts. All study data is stored in a database for easy exporting to statistical packages. Through the use of the PAS Reporter utility, zone entries, time in zone, distance traveled, speed and resting time can also be reported. Powerful software features provide the ability to graphically view the subject's path.

• **Maximum Stations: 16 with or without rearing frames and floor pokes** (For older PCI interface PAS Open Field Systems please contact San Diego Instruments to determine if your system can be supported, expanded and/or converted to a USB interface.)

2325-0409	PHOTOBREAM ACTIVITY SYSTEM - OPEN FIELD Includes: Control Unit with USB interface, software, 16 x 16 photobeam mounting frame and all cables & connectors.
7500-0334	Photobeam Activity System - Open Field Additional Test Station
2325-0410	PHOTOBREAM ACTIVITY SYSTEM - OPEN FIELD - 4 HOLE FLOOR POKE SYSTEM Includes: Control Unit with USB interface, software, 16 x 16 photobeam mounting frame, 4 hole floor poke and all cables & connectors
7500-0336	Photobeam Activity System - Open Field - 4 Hole Poke Additional Test Station
2325-0411	PHOTOBREAM ACTIVITY SYSTEM - OPEN FIELD - 16 HOLE FLOOR POKE SYSTEM Includes: Control Unit with USB interface, software, 16 x 16 photobeam mounting frame, 16 hole floor poke and all cables & connectors.
7500-0337	Photobeam Activity System - Open Field - 16 Hole Poke Additional Test Station


PHOTOBREAM ACTIVITY SYSTEM - OPEN FIELD OPTIONAL COMPONENTS & ACCESSORIES

7500-0335	Photobeam Activity System - Open Field Rearing Frame
7001-0364	Photobeam Activity System - Dark Box Insert
6300-0100	Photobeam Activity System - Open Field Windows® Software Version Upgrade
6300-0107	Photobeam Activity System - Open Field FDA Part 11 Compliance Software
6603-0101	Photobeam Activity System - Open Field GLP Validation Document Set
2325-0104	Photobeam Activity System - Open Field Isolation Cabinet

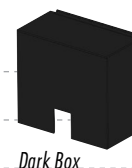

Dark Box

ROTOMETER SYSTEM

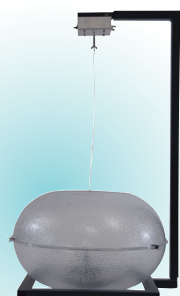

*Optional Rotometer
System Mounting Frame*

Rotometer Test Station

SDI's Rotometer System is powerful activity system designed to record the subject's clockwise and counter clockwise rotations. SDI's Rotometer has been used widely in Traumatic Brain Injury and Spinal Cord Injury studies. Each station starts and runs independently making continuous testing quick and easy. An optional quarter turn version of the system can be utilized to count partial rotations. SDI's Rotometer takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use.


• **Maximum Stations: 16**

(For older PCI interface Rotometer Systems please contact San Diego Instruments to determine if your system can be supported, enlarged and/or converted to a USB interface.)


2325-0412	ROTOMETER SYSTEM Includes: Control Unit with USB interface, software, test station and all cables & connectors.
7500-0338	Rotometer System - Additional Test Station
2325-0413	ROTOMETER SYSTEM WITH QUARTER TURN ADD-ON FEATURE Includes: Control Unit with USB interface, software, test station and all cables & connectors.
7500-0339	Rotometer System - Additional Test Station with Quarter Turn Add-on Feature
	ROTOMETER SYSTEM OPTIONAL COMPONENTS & ACCESSORIES
7500-0279	Rotometer Test Station Mounting Frame
7500-0304	Rotometer Rat Enclosure
7500-0305	Rotometer Mouse Enclosure
6300-0300	Rotometer Software Version Upgrade

PLACE PREFERENCE SYSTEM

PLACE PREFERENCE

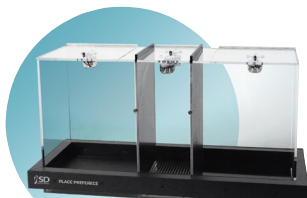

Place Preference 3-Compartment System
with Photobeams

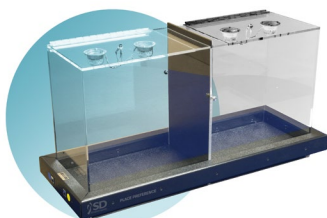

Place Preference 2-Compartment System
with Photobeams

SDI's Place Preference System is the most powerful, flexible and easy to use system for measuring place preference. SDI's Place Preference System employs a 4 x 16 photobeam array to precisely record entries into chamber and time in chamber. SDI's Place Preference System takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. All study data is stored in a database for easy exporting to statistical packages. SDI's Place Preference System can also be utilized as an activity system. The Place Preference Reporter reduces the collected data to show the number of entries in a compartment and the time spent for each entry.

SDI's Place Preference test enclosure is designed with clear acrylic walls providing the ability to attach any type of cue. Removable floors allow the user to create custom floors from virtually any material to meet specific study requirements. High intensity LED ceiling lights are long lasting and emit no heat. Manually operated doors rise fully, providing an unobstructed pathway for the subject to cross between chambers.

• Maximum Stations: 16

(For older PCI interface Place Preference Systems please contact San Diego Instruments to determine if your system can be supported, expanded and/or converted to a USB interface.)

PLACE PREFERENCE SYSTEM

Includes: Control Unit with USB interface, software, 4 x 16 photobeam configuration, LED lights, removable floors and all cables & connectors.

2325-0414	Place Preference System (3 Compartment)
7500-0340	Place Preference System - Additional Test Station (3 Compartment)
2325-0415	Place Preference System (2 Compartment)
7500-0341	Place Preference System - Additional Test Station(2 Compartment)
2325-0416	Place Preference System (Combo 2 or 3 Compartment)
7500-0342	Place Preference System - Additional Test Station(Combo 2 or 3 Compartment)

PLACE PREFERENCE ENCLOSURES FOR USE WITH HUMAN OBSERVATION OR VIDEO TRACKING

(No photobeams and no LED lights in the tops)

7500-0241	Place Preference Enclosure (3 Compartment)
7500-0293	Place Preference Enclosure (2 Compartment)
7500-0294	Place Preference Enclosure (Combo 2 or 3 Compartment)

TREMOR MONITOR™ SYSTEM

TREMOR MONITOR™



Tremor Monitor Station

The Tremor Monitor™ System provides the ability to record, review, analyze and report rapid repetitive movements in rats and mice. Tremor Monitor takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. Tremor Monitor provides the ability to configure multiple test stations, allowing rapid testing of a large number of subjects. Tremor Monitor accurately differentiates tremor events from ambulatory/stereotyped movements. Tremor Monitor utilizes an ultra-sensitive movement sensor to record continuous movement waveforms at 128Hz for up to 30 minutes. The user can review the results in adjustable interval lengths.

• **Maximum Stations: 8** (For older PCI interface Tremor Monitor Systems please contact San Diego Instruments to determine if your system can be supported, expanded and/or converted to a USB interface.)

2325-0424

TREMOR MONITOR SYSTEM

Includes: 1 ABS cabinet, 1 animal enclosure, software, connection chassis with USB interface and all cables & connectors.

2325-0425

Tremor Monitor System - Additional Test Station

Includes: 1 ABS cabinet, 1 animal enclosure and all cables & connectors.

2325-0246

TREMOR MONITOR/SR-LAB ABS COMBO SYSTEM

Includes: 1 ABS cabinet, 1 Tremor Monitor animal enclosure, 1 SR-LAB animal enclosure, Tremor Monitor software, SR-LAB software, connection chassis with USB interface and all cables & connectors.

2325-0427

Tremor Monitor/SR-LAB ABS Combo System - Additional Test Station

Includes: 1 ABS cabinet, 1 Tremor Monitor animal enclosure and 1 SR-LAB animal enclosure and all cables & connectors.

7500-0348


Tremor Monitor Animal Enclosure

2325-0287

Tremor Standardization Unit

ANY-MAZE™ SYSTEM

ANY-MAZE™

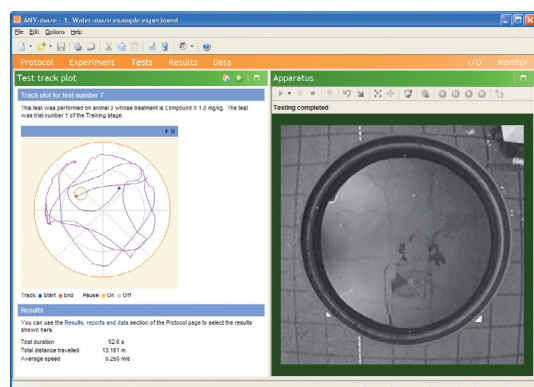

ANY-maze™ is a flexible video tracking system designed to automate testing in behavioral experiments. Packed with advanced features ANY-maze is one of the most comprehensive video tracking systems available today. ANY-maze can perform tests in up to sixteen pieces of apparatus simultaneously. ANY-maze provides multiple measurements in these key areas:

- Apparatus measures
- Key (lever) measures
- Zone measures
- Input switch measures
- Sequence measures

The standard ANY-maze bundle includes software, camera, camera ceiling mount and remote switch. Study specific bundles are also available that add unique components to the standard bundle. When contacting San Diego Instruments for a bundle including an apparatus please include whether you will be testing rats, mice or both and your desired color for the apparatus (see the Mazes and Enclosures section for the color choices for each apparatus).

ANY-MAZE™ STANDARD BUNDLES

2325-0275	ANY-maze Video Tracking System Bundle
2325-8275	ANY-maze Video Tracking System Bundle with Laptop
2325-9275	ANY-maze Video Tracking System Bundle with Cobalt Desktop

ANY-MAZE™ STANDARD BUNDLES WITH APPARATUS

2325-7235	ANY-maze Video Tracking System Bundle with Mouse Barnes Maze
2325-7236	ANY-maze Video Tracking System Bundle with Rat Barnes Maze
2325-7237	ANY-maze Video Tracking System Bundle with Mouse Y-Maze
2325-7238	ANY-maze Video Tracking System Bundle with Rat Y-Maze
2325-7239	ANY-maze Video Tracking System Bundle with Mouse T-Maze
2325-7240	ANY-maze Video Tracking System Bundle with Rat T-Maze
2325-7241	ANY-maze Video Tracking System Bundle with Mouse Elevated Plus-Maze
2325-7242	ANY-maze Video Tracking System Bundle with Rat Elevated Plus-Maze
2325-7243	ANY-maze Video Tracking System Bundle with Mouse Zero Maze
2325-7244	ANY-maze Video Tracking System Bundle with Rat Zero Maze
2325-7245	ANY-maze Video Tracking System Bundle with Mouse Radial Arm Maze
2325-7246	ANY-maze Video Tracking System Bundle with Rat Radial Arm Maze
2325-7247	ANY-maze Video Tracking System Bundle with Mouse Water Maze Components <i>See note #2</i>
2325-7248	ANY-maze Video Tracking System Bundle with Rat Water Maze Components <i>See note #2</i>
2325-7249	ANY-maze Video Tracking System Bundle with Place Preference
2325-7250	ANY-maze Video Tracking System Bundle with Multiple Unit Open Field
2325-7251	ANY-maze Video Tracking System Bundle with Single Unit Open Field
2325-7252	ANY-maze Video Tracking System Bundle with Forced Swim Test Enclosure
2325-7255	ANY-maze Video Tracking System Bundle with Watertight Mouse Radial Arm Maze
2325-7256	ANY-maze Video Tracking System Bundle with Watertight Rat Radial Arm Maze
2325-8235	ANY-maze Video Tracking System Bundle with Laptop and Mouse Barnes Maze
2325-8236	ANY-maze Video Tracking System Bundle with Laptop and Rat Barnes Maze
2325-8237	ANY-maze Video Tracking System Bundle with Laptop and Mouse Y-Maze
2325-8238	ANY-maze Video Tracking System Bundle with Laptop and Rat Y-Maze
2325-8239	ANY-maze Video Tracking System Bundle with Laptop and Mouse T-Maze
2325-8240	ANY-maze Video Tracking System Bundle with Laptop and Rat T-Maze
2325-8241	ANY-maze Video Tracking System Bundle with Laptop and Mouse Elevated Plus Maze
2325-8242	ANY-maze Video Tracking System Bundle with Laptop and Rat Elevated Plus Maze

STANDARD BUNDLES WITH APPARATUS (CONT.)

2325-8243	ANY-maze Video Tracking System Bundle with Laptop and Mouse Zero Maze
2325-8244	ANY-maze Video Tracking System Bundle with Laptop and Rat Zero Maze
2325-8245	ANY-maze Video Tracking System Bundle with Laptop and Mouse Radial Arm Maze
2325-8246	ANY-maze Video Tracking System Bundle with Laptop and Rat Radial Arm Maze
2325-8247	ANY-maze Video Tracking System Bundle with Laptop and Mouse Water Maze Components <i>See note #2</i>
2325-8248	ANY-maze Video Tracking System Bundle with Laptop and Rat Water Maze Components <i>See note #2</i>
2325-8249	ANY-maze Video Tracking System Bundle with Laptop and Place Preference
2325-8250	ANY-maze Video Tracking System Bundle with Laptop and Multiple Unit Open Field
2325-8251	ANY-maze Video Tracking System Bundle with Laptop and Single Unit Open Field
2325-8252	ANY-maze Video Tracking System Bundle with Laptop and Forced Swim Test Enclosure
2325-8253	ANY-maze Video Tracking System Bundle with Laptop and Submersible Mouse Radial Arm Maze
2325-8254	ANY-maze Video Tracking System Bundle with Laptop and Submersible Rat Radial Arm Maze
2325-8255	ANY-maze Video Tracking System Bundle with Laptop and Watertight Mouse Radial Arm Maze
2325-8256	ANY-maze Video Tracking System Bundle with Laptop and Watertight Rat Radial Arm Maze
2325-9235	ANY-maze Video Tracking System Bundle with Desktop and Mouse Barnes Maze
2325-9236	ANY-maze Video Tracking System Bundle with Desktop and Rat Barnes Maze
2325-9237	ANY-maze Video Tracking System Bundle with Desktop and Mouse Y-Maze
2325-9238	ANY-maze Video Tracking System Bundle with Desktop and Rat Y-Maze
2325-9239	ANY-maze Video Tracking System Bundle with Desktop and Mouse T-Maze
2325-9240	ANY-maze Video Tracking System Bundle with Desktop and Rat T-Maze
2325-9241	ANY-maze Video Tracking System Bundle with Desktop and Mouse Elevated Plus-Maze
2325-9242	ANY-maze Video Tracking System Bundle with Desktop and Rat Elevated Plus-Maze
2325-9243	ANY-maze Video Tracking System Bundle with Desktop and Mouse Zero Maze
2325-9244	ANY-maze Video Tracking System Bundle with Desktop and Rat Zero Maze
2325-9245	ANY-maze Video Tracking System Bundle with Desktop and Mouse Radial Arm Maze
2325-9246	ANY-maze Video Tracking System Bundle with Desktop and Rat Radial Arm Maze
2325-9247	ANY-maze Video Tracking System Bundle with Desktop and Mouse Water Maze Components <i>See note #2</i>
2325-9248	ANY-maze Video Tracking System Bundle with Desktop and Rat Water Maze Components <i>See note #2</i>
2325-9249	ANY-maze Video Tracking System Bundle with Desktop and Place Preference
2325-9250	ANY-maze Video Tracking System Bundle with Desktop and Multi-Unit Open Field
2325-9251	ANY-maze Video Tracking System Bundle with Desktop and Single Unit Open Field
2325-9252	ANY-maze Video Tracking System Bundle with Desktop and Forced Swim Test Enclosure

2325-9255	ANY-maze Video Tracking System Bundle with Desktop and Watertight Mouse Radial Arm Maze
2325-9256	ANY-maze Video Tracking System Bundle with Desktop and Watertight Rat Radial Arm Maze
ANY-MAZE SYSTEM OPTIONAL COMPONENTS & ACCESSORIES	
7500-0290	Infrared Camera with Illuminator for Desktop
7500-0295	Infrared Camera with Illuminator for Laptop
2500-0514	Camera Ceiling Mount
5100-0058	Remote Switch
7500-0311	Video Camera Bundle(Includes: Video Camera(USB), Varifocal Lens, Ceiling Mount and 16 ft(5 M) USB Cable)

Note #2: Includes Water tank, Escape Platform, and Drain Kit

MAZES & ENCLOSURES

MAZES & ENCLOSURES



SDI's Mazes and Enclosures are made of durable, high density, non-porous plastic. Sizes for rats and mice are available. Unless otherwise noted, they are available in white, black and beige. SDI Mazes and Enclosures are optimally designed for use with SDI ANY-Maze™ Video Tracking System. (see pages 14-15)

Water Maze Tank, Elevated Plus Maze, Single Unit Open Field Enclosure

ZERO MAZE

2325-0230	Zero Maze - Black (Mouse)
2325-0229	Zero Maze - White (Mouse)
2325-0231	Zero Maze - Beige (Mouse)
2325-0227	Zero Maze - Black (Rat)
2325-0226	Zero Maze - White (Rat)
2325-0228	Zero Maze - Beige (Rat)

ELEVATED PLUS MAZE

7001-0225	Elevated Plus Maze - Black (Mouse)
7001-0336	Elevated Plus Maze - White (Mouse)
7001-0316	Elevated Plus Maze - Beige (Mouse)
7001-0227	Elevated Plus Maze - Black (Rat)
7001-0338	Elevated Plus Maze - White (Rat)
7001-0337	Elevated Plus Maze - Beige (Rat)

Y-MAZE

Includes: Interchangeable floors (plain and pellet holes), manual doors at each upper arm and starter door in base leg

7001-0422	Y-Maze - Black - Mouse
7001-0419	Y-Maze - White - Mouse
7001-0425	Y-Maze - Beige - Mouse
7001-0423	Y-Maze - Clear Walls - Black Floors - Mouse
7001-0420	Y-Maze - Clear Walls - White Floors - Mouse
7001-0424	Y-Maze - Clear Walls - Beige Floors - Mouse
7001-0426	Y-Maze - Clear Walls - Clear Floors - Mouse
7001-0431	Y-Maze - Black - Rat
7001-0428	Y-Maze - White - Rat
7001-0434	Y-Maze - Beige - Rat
7001-0432	Y-Maze - Clear Walls - Black Floors - Rat
7001-0429	Y-Maze - Clear Walls - White Floors - Rat
7001-0433	Y-Maze - Clear Walls - Beige Floors - Rat
7001-0435	Y-Maze - Clear Walls - Clear Floors - Rat

T-MAZE

Includes: Interchangeable floors(plain and pellet holes) and manual doors at each upper arm and starter door in base leg

7001-0404	T-Maze - Black - Mouse
7001-0401	T-Maze - White - Mouse
7001-0407	T-Maze - Beige - Mouse
7001-0405	T-Maze - Clear Walls - Black Floors - Mouse
7001-0402	T-Maze - Clear Walls - White Floors - Mouse
7001-0406	T-Maze - Clear Walls - Beige Floors - Mouse
7001-0408	T-Maze - Clear Walls - Clear Floors - Mouse
7001-0413	T-Maze - Black - Rat
7001-0410	T-Maze - White - Rat
7001-0416	T-Maze - Beige - Rat
7001-0414	T-Maze - Clear Walls - Black Floors - Rat
7001-0411	T-Maze - Clear Walls - White Floors - Rat
7001-0415	T-Maze - Clear Walls - Beige Floors - Rat
7001-0417	T-Maze - Clear Walls - Clear Floors - Rat

RADIAL ARM MAZE

Includes: Interchangeable floors (plain and pellet holes) and manual doors at each arm

7001-0440	Radial Arm Maze - Black - Mouse
7001-0437	Radial Arm Maze - White - Mouse
7001-0443	Radial Arm Maze - Beige - Mouse
7001-0441	Radial Arm Maze - Clear Walls Black Floors - Mouse
7001-0438	Radial Arm Maze - Clear Walls White Floors - Mouse
7001-0442	Radial Arm Maze - Clear Walls Beige Floors - Mouse
7001-0444	Radial Arm Maze - Clear Walls Clear Floors - Mouse
7001-0449	Radial Arm Maze - Black - Rat
7001-0446	Radial Arm Maze - White - Rat
7001-0452	Radial Arm Maze - Beige - Rat
7001-0450	Radial Arm Maze - Clear Walls Black Floors - Rat
7001-0447	Radial Arm Maze - Clear Walls White Floors - Rat
7001-0451	Radial Arm Maze - Clear Walls Beige Floors - Rat
7001-0453	Radial Arm Maze - Clear Walls Clear Floors - Rat

BARNES MAZE

7001-0235	Barnes Maze 36" - White (Mouse)
7001-0324	Barnes Maze 36" - Black (Mouse)
7001-0340	Barnes Maze 36" - Beige (Mouse)
7001-0326	Barnes Maze 48" - White (Rat)
7001-0325	Barnes Maze 48" - Black (Rat)
7001-0341	Barnes Maze 48" - Beige (Rat)

MULTIPLE UNIT OPEN FIELD

7001-0353	Multiple Unit Open Field - Black
7001-0354	Multiple Unit Open Field - White
7001-0067	Multiple Unit Open Field - Beige

SINGLE UNIT OPEN FIELD


7001-0355	Single Unit Open Field - Black
7001-0356	Single Unit Open Field - White
7001-0357	Single Unit Open Field - Beige

WATER TEST COMPONENTS


SDI's Morris Water Maze components include a tank, escape platform and drain kit. Sizes for rats and mice are available. The water tank is made of durable, high density polyethylene and comes in either white or blue. The Escape Platform is available with either a round or square multi-step top in clear acrylic. The Escape Platform includes a dowel rod with colored tape that can be inserted in the top to provide a flag for training. The optional drain kit eliminates the need for drain piping. SDI's Morris Water Maze components are optimally designed for use with SDI ANY-Maze™ Video Tracking System. (see pages 14-15)

Morris

7500-0271	Water Maze Platform - Clear Square Top
7500-0274	Water Maze Platform - Clear Round Top
7500-0309	Water Maze Platform Adjustable Height (11" to 22") - Clear Square Top
7500-0310	Water Maze Platform Adjustable Height (11" to 22") - Clear Round Top
7000-0718	Water Rat Tank - 6 ft. - Blue
7000-0719	Water Rat Tank - 6 ft - White
7000-0720	Water Mouse Tank - 4 ft - Blue
7000-0723	Water Mouse Tank - 4 ft - White
7500-0217	Water Tank Drain Kit

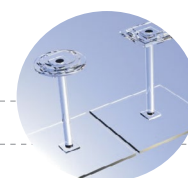

Water Maze Platforms
Fixed

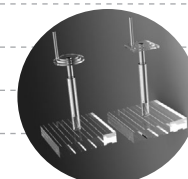


Water Maze Platforms
Adjustable



Bulkhead fitting (pre-mounted on tank) provides drain hole with 3/4" FIPT to connect Drain Kit or user supplied drain piping

Water Test Radial Arm Mazes

7001-0461	Radial Arm Maze - Watertight - Mouse - Black
7001-0462	Radial Arm Maze - Watertight - Mouse - White
7001-0463	Radial Arm Maze - Watertight - Mouse - Beige -
7001-0464	Radial Arm Maze - Watertight - Rat - Black
7001-0465	Radial Arm Maze - Watertight - Rat - White
7001-0466	Radial Arm Maze - Watertight - Rat - Beige
7001-0398	Radial Arm Maze - Watertight Escape Platform - Mouse
7001-0399	Radial Arm Maze - Watertight Escape Platform - Rat

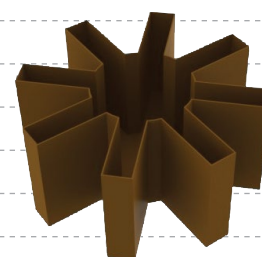

Radial Arm Maze- Watertight

HUMAN TESTING SYSTEMS

SR-HLAB™

SR-HLAB™

Eyeblink Conditioning System

ERP-LAB™

Human Testing Combination Systems



SR-HLAB EMG System

The SR-HLAB™ Startle Response System for humans measures startle response by recording the eyeblink reflex via electromyographic (EMG) or photoelectric cell (PEC) recording. SR-HLAB is a complete hardware and software solution for a wide variety of startle applications. Intuitive yet powerful features cater to users from the educated non-specialist to the sophisticated behaviorist. SR-HLAB is portable, making it the first startle analysis system that can be easily transported from the laboratory to offsite test facilities. The complete system weighs only 8 pounds and comes with a carrying case.

SR-HLAB takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use. SR-HLAB is available in two models: SR-HLAB EMG and SR-HLAB PEC. The EMG model uses electrodes positioned over the subject's Orbicularis Oculi muscle to capture the eyeblink reflex. The PEC model uses a photoelectric cell monitoring the eyelid to capture the eyeblink reflex. Both models utilize the same software and scoring algorithms.

2325-0417 **SR-HLAB SYSTEM EMG** (Electromyographic sensing for eyeblink)

2325-0428 **SR-HLAB SYSTEM PEC** (Photoelectric Cell sensing for eyeblink)

OPTIONAL COMPONENTS & ACCESSORIES

7500-0343 EMG Response Module

7500-0344 PEC Response Module (for photo electric cell sensor or heart rate sensor)

7500-0345 Galvanic Skin Resistance Module

7500-0346 Digital I/O Module

2325-0119 Tactile Kit Option – air puff with user supplied air tank

2325-0253 Portable Airpuff Unit – small compressor supplies air

2200-0092 Relay Chassis (4 - 12 volt Solid State Relays)

7500-0198 Slave Mouse

6200-0189 Heart Rate Sensor Unit

6200-0092 EMG Grey Yoke 1 Channel Cable

1900-0008 EMG Reusable Electrodes (3-pack with red/white/black lead wires)

1900-0005 EMG Reusable Electrodes-Collars (5 packs of 100 each, min 500)

1900-0006 EMG Reusable Electrodes-Saline Gel (Minimum Purchase 5 tubes)

1900-0009 EMG Disposable Electrodes (3 per pack, min 20 packs)

1900-0011 EMG Disposable Electrodes-Lead Wire Set (red/white/black)

1900-0010 EMG Skin Preparation Abrasive Pads

6300-0900 SR-HLAB Windows Software Version Upgrade

EYEBLINK CONDITIONING SYSTEM™



Eyeblink Conditioning System

SDI's Eyeblink Conditioning System is used widely in basic learning studies, psychopharmacology research and the examination of the behavioral consequences of Alzheimer's and other diseases. SDI's Eyeblink Conditioning System employs an integrated photoelectric cell with air tube to monitor the eyeblink reflex and deliver the air puff under software control. For complete portability the Portable Air Puff Unit with built in compressor is available, eliminating the need for an air tank. SDI's Eyeblink Conditioning System takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use.

2325-0418

EYEBLINK CONDITIONING SYSTEM

Includes: Tactile Kit Option – air puff with user supplied air tank

EYEBLINK OPTIONAL COMPONENTS & ACCESSORIES

2325-0253	Portable Airpuff Unit –small compressor supplies air
2325-0119	Tactile Kit Option – air puff with user supplied air tank
7500-0209	Headset - Fleece Band Infrared Sensor & Air Tube
7500-0347	Headset - Fleece Band Infrared Sensor
7500-0213	Headset - Fleece Bands Only (3 pack)
7500-0220	Visual Indicator Unit
7500-0343	EMG Response Module
7500-0344	PEC Response Module(for photo electric cell sensor or heart rate sensor)
7500-0345	Galvanic Skin Resistance Module
7500-0346	Digital I/O Module
2200-0092	Relay Chassis (4 - 12 volt Solid State Relays)
7500-0198	Slave Mouse
6200-0189	Heart Rate Sensor Unit
6200-0092	EMG Grey Yoke 1 Channel Cable
1900-0008	EMG Reusable Electrodes(3-pack with red/white/black lead wires)
1900-0005	EMG Reusable Electrodes-Collars (5 packs of 100 each, min 500)
1900-0006	EMG Reusable Electrodes-Saline Gel (Minimum Purchase 5 tubes)
1900-0009	EMG Disposable Electrodes(3 per pack, min 20 packs)
1900-0011	EMG Disposable Electrodes-Lead Wire Set(red/white/black)
1900-0010	EMG Skin Preparation Abrasive Pads
6300-1000	Eyeblink Conditioning System Windows Software Version Upgrade

ERP-LAB™



ERP-LAB System

The ERP-LAB™ System is a scaled down or “light” version of the traditional Evoked Response Potential (ERP) test system offering P50, P300 or Mismatch Negativity. ERP-LAB is designed for neuropsychiatric researchers who are specifically interested in testing already identified landmark areas of the brain. ERP-LAB utilizes two channels to assure accurate study data. The EEG channel records the P50, P300 or Mismatch Negativity responses while the EOG channel is used for artifact detection. ERP-LAB takes full advantage of the Windows® operating system with data organization and management software that combines power and flexibility with ease of use.

2325-0419

ERP-LAB SYSTEM

ERP-LAB OPTIONAL COMPONENTS & ACCESSORIES

7500-0344	PEC Response Module (for photo electric cell sensor or heart rate sensor)
7500-0345	Galvanic Skin Resistance Module
7500-0346	Digital I/O Module
6200-0189	Heart Rate Sensor Unit
6200-0092	EMG Grey Yoke 1 Channel Cable
6200-0093	EMG Reusable Electrodes(3-pack with red/white/black lead wires)
1900-0005	EMG Reusable Electrodes-Collars (5 packs of 100 each, min 500)
1900-0006	EMG Reusable Electrodes-Saline Gel (Minimum Purchase 5 tubes)
1900-0009	EMG Disposable Electrode (3 per pack, min 20 packs)
1900-0011	EMG Disposable Electrodes-Lead Wire Set (red/white/black)
1900-0010	EMG Skin Preparation Abrasive Pads
6300-1100	ERP-LAB Windows Software Version Upgrade

HUMAN TESTING COMBINATION SYSTEMS

The three human testing systems SR-HLAB, Eyeblink Conditioning and ERP-LAB use the same hardware thus they can be offered as combination systems. In a combination system all software is included to allow the user to switch between test types. All the requisite response and stim modules are included. Other modules are included as required based on the selection of an Optional Component or Accessory selected from the lists following the individual systems.

2325-0420

SR-HLAB SYSTEM/EYEBLINK CONDITIONING SYSTEM

Includes: Stim Module, EMG Response Module, PEC Response Module, EMG Electrodes, Fleece Head Band with PEC sensor, Tactile Kit, SR-HLAB software and Eyeblink Conditioning Software

2325-0421

SR-HLAB SYSTEM/ERP-LAB SYSTEM

Includes: Stim Module, 2-EMG Response Module, 2 sets of EMG Electrodes, SR-HLAB software and ERP-LAB Software.

2325-0422

EYEBLINK CONDITIONING SYSTEM/ERP-LAB SYSTEM

Includes: Stim Module, 2-EMG Response Module, 1-PEC Response Module, 2 sets of EMG Electrodes, Fleece Head Band with PEC sensor, Tactile Kit, Eyeblink Conditioning Software and ERP-LAB software

2325-0423

SR-HLAB SYSTEM/EYEBLINK CONDITIONING SYSTEM/ERP-LAB SYSTEM

Includes: Stim Module, 2-EMG Response Module, 1-PEC Response Module, 2 sets of EMG Electrodes, Fleece Head Band with PEC sensor, Tactile Kit, SR-HLAB Software, Eyeblink Conditioning Software and ERP-LAB software

TEST ACCESSORIES

TEST ACCESSORIES



AIRSTIM™

AIRSTIM™ with Software

Portable Air Puff Unit

Portable Air Puff with Software

Animal Stimulator

Animal Shocker

Shock Level Tester

White Noise Generator

2325-0129

AIRSTIM UNIT

7500-0207

AIRSTIM AIR DELIVERY TUBES

2325-0130

AIRSTIM UNIT WITH SOFTWARE CONTROL

2325-0370

PORTABLE AIR PUFF UNIT

2325-0371

PORTABLE AIR PUFF UNIT WITH SOFTWARE CONTROL

2325-0360

ANIMAL STIMULATOR (shock, light and sound under software control)

2325-0124

PROGRAMMABLE ANIMAL SHOCKER

2325-0358

SHOCK LEVEL TESTER

2325-0144

WHITE NOISE GENERATOR

White Noise Generator for Lab Use

Includes: 24 VDC power supply and 1 speaker (volume up to 88dB)



Animal Stimulator



Programmable Animal Shocker



White Noise Generator

SDI GENERAL INFO

TO REQUEST A QUOTE OR PLACE AN ORDER

San Diego Instruments offers several convenient ways to get a quote or order:

- To get a quote online: sandiegoinstruments.com/quote
- To get a quote or order by email: sales@sandiegoinstruments.com
- To get a quote or order by phone: 858-530-2600
- To get a quote or order by fax: 858-530-2646

TO SPEAK WITH A SPECIALIST

Our pre-sales consultation process is designed to ensure that you purchase the right system for your unique research demands. The choice among system configurations can be a complex one. We want to hear about your study so that we may best assist you with the decision. If you have any questions or would like to speak to one of our specialists about your specific needs, please call: 858-530-2600

TO LEARN MORE ABOUT SDI BEHAVIORAL TESTING SYSTEMS

Our website offers a complete source of product information including features and benefits, technical specifications, data sheets, scholarly citations and product overview presentations. To access and download this information, please visit: www.sandiegoinstruments.com

TO VIEW PRODUCT OVERVIEW PRESENTATIONS

Our product overview presentations provide an engaging way to discover the wide range of features and benefits provided by SDI systems. These short (5-10 minutes) narrated presentations provide a comprehensive insight into the considerable advantages of SDI systems. Each presentation covers topics such as setting up an experiment, diagnostics routines, running a session, viewing reports and exporting data. To view product presentations online, please visit: sandiegoinstruments.com/presentations



Designed by Scientists for Scientists

LAPTOP COMPATIBLE SYSTEMS

SR-LAB™, GEMINI™, ROTOR-ROD™, Tremor Monitor™, Freeze Monitor™, Grip Strength, Tail Flick, Hotplate, Hot/Cold Plate Pressure Analgesia Meter, SR-HLAB™, ERP-LAB™ and Eyeblink Conditioning System.

MINIMUM COMPUTER REQUIREMENTS

Windows XP/Windows 7 computer system with USB connection and PCI slot(where noted). Minimum disk and memory size specified to support Windows XP/Windows 7 are acceptable.

TERMS OF SALE

All San Diego Instruments' products are subject to software, firmware and design copyrights and trademarks. Please refer to you User's Manual for restrictions on duplication and other usage of San Diego Instruments' products.

WARRANTY

One-year factory warranty on all parts and labor. Missing parts must be reported within 30 days of receipt to be covered under warranty. Please refer to San Diego Instruments' Warranty Statement for details.

DISCLAIMER

The equipment described in this catalog is designed for research and educational applications and is not intended for the diagnosis, alleviation, treatment, monitoring or prevention of disease, injury or handicap.



San Diego Instruments, Inc.
9155 Brown Deer Rd, Suite 8
San Diego, CA 92121
Ph: 858-530-2600
Fax: 858-530-2646
www.sandiegoinstruments.com

©2014 San Diego Instruments. All rights reserved. Prices and specifications may be changed without notification.



Mixed Sources
Product group from well-managed
forests and other controlled sources
www.fsc.org Cert no. SW-COC-001530
© 1996 Forest Stewardship Council